

M38 Mobile Computing Models



Non-Powered Computer Cart



Powered Computer Cart



Medication Expansion Pack

Mobile Computing Carts

NON-POWERED CARTS

- Configurations starting at 70 lbs
- Accommodates up to 17" notebook and tablet PC's
- External spiral cord recharges notebook
- Will accommodate LCD monitors if used in applications where cart is always plugged in

ALL POWERED CARTS

- Batteries: SLA AGM 35 or 55 AMP or Lithium Iron Phosphate (LiFePO4) 40 AMP*
- Certified to UL 60601-1 standard for patient point of contact
- Input 120V 60 Hz; Charger: 10 AMP SLA and 20 AMP LiFePO4
- Automatic transfer switch enables full use of PC without affecting battery charge time
- Gentle shut-down at critically low power levels
- Power indicator at work-surface height with audible alarm
- Accepts LCD monitors up to 20"
- Auto shut down of PC prevents deep discharge
- Software tracks battery age and notifies when replacements are due
- Email setup sends cart specific critical alerts

AC POWERED CARTS

- Configurations starting at 104 lbs
- External CPU configuration accommodates 2 ¾" to 4" wide CPU's

DC POWERED CARTS


- High efficiency system provides more than 12 hours of run-time**
- Configurations starting at 103 lbs
- 3 DC adapter power outlets
- Supplemental on-screen battery indicates time remaining

LITHIUM IRON PHOSPHATE (LiFePO4) BATTERY

- 40 AMP LiFePO4 provides 5 times longer life than 55 AH SLA
- Charging Time: 2 Hours
- Weight: 14 lbs (65% lighter than SLA equivalent)

MEDICATION EXPANSION PACK

- Secure drawer system can be ordered factory installed or added at a future date
- Electric Keyless Entry PIN Code System: 0-9 digits, 128 codes per cart, update and audit PIN codes remotely
- Drawer lock time-out (seconds): 1-255
- Drawer system is reconfigurable on-the-go

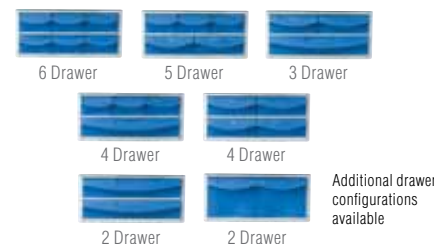

*LiFePO4 for DC Carts available Q2 2011
 **Run times will vary depending on technology used and charging patterns

All M38 Models

GENERAL

- Base size: 17" x 17" contoured to 14.5" wide at center
- Height adjustment: 16" adjustable tension spring
- Contoured front and rear handle
- Work Surface size: 21" w x 16" d
- LED Keyboard light with auto shut-off
- Accommodates 8 ¾" d x 18" w keyboards
- Keyboard tray adjustment: 6" height; +10° to -15° angle; 180° swivel
- Keyboard tray has left/right mouse pads, wrist rest, and easy access mouse holder
- Work surface Document Protector: 16" w x 9.5" d
- Bin Accommodates barcode scanner heads up to 4 x 3 ¾"
- Casters: 4" precision bearing casters/2 locking
- T-track provides universal mounting point for accessories
- Smooth, cleanable surfaces
- High touch surfaces made with antimicrobial materials
- 2.5 ft hospital grade spiral cord extends to 8 ft
- Internal cable management
- Internal CPU storage: 15 ¼" w x 13 ¼" d x 3 ¼" h
- 3 available power outlets

STORAGE

- 2 general purpose storage bins
- Storage for cables under work surface
- Drawer with mechanical keyless lock (optional)
- Wire basket (optional)

WARRANTY

- 3 year mechanical; 2 year electrical
- Extended warranty available (please inquire)

OPTIONS

- 35 or 55 AMP SLA battery or 40 AMP LiFePO4 battery
- Wire basket, side bins, IV Poles, vital signs mounts, scanner/printer shelf
- Drawer with mechanical keyless PIN lock
- Universal mounting bracket to add items such as bins, sharps containers, hand sanitizers, and IV poles

SUITED FOR

- Bedside documentation
- Clinical documentation
- eMAR/Bar-coding implementations
- Automated dispensing cabinets such as Pyxis and Accudose
- Patient registration
- Physician order entry/rounds
- Med/Surg, Emergency, Ambulatory/Outpatient, OR's, ICU, Clinics



SERVICE GUARANTEE

We are committed to providing best in class service to maximize your uptime. Our service guarantee includes:


- 2 hour response time
- Parts shipped next business day
- Help desk support
- Dedicated Technical Account Manager



Visit Rubbermaidmedical.com or call 888-859-8294 to see our complete line of medical products



Mobile Computing Carts

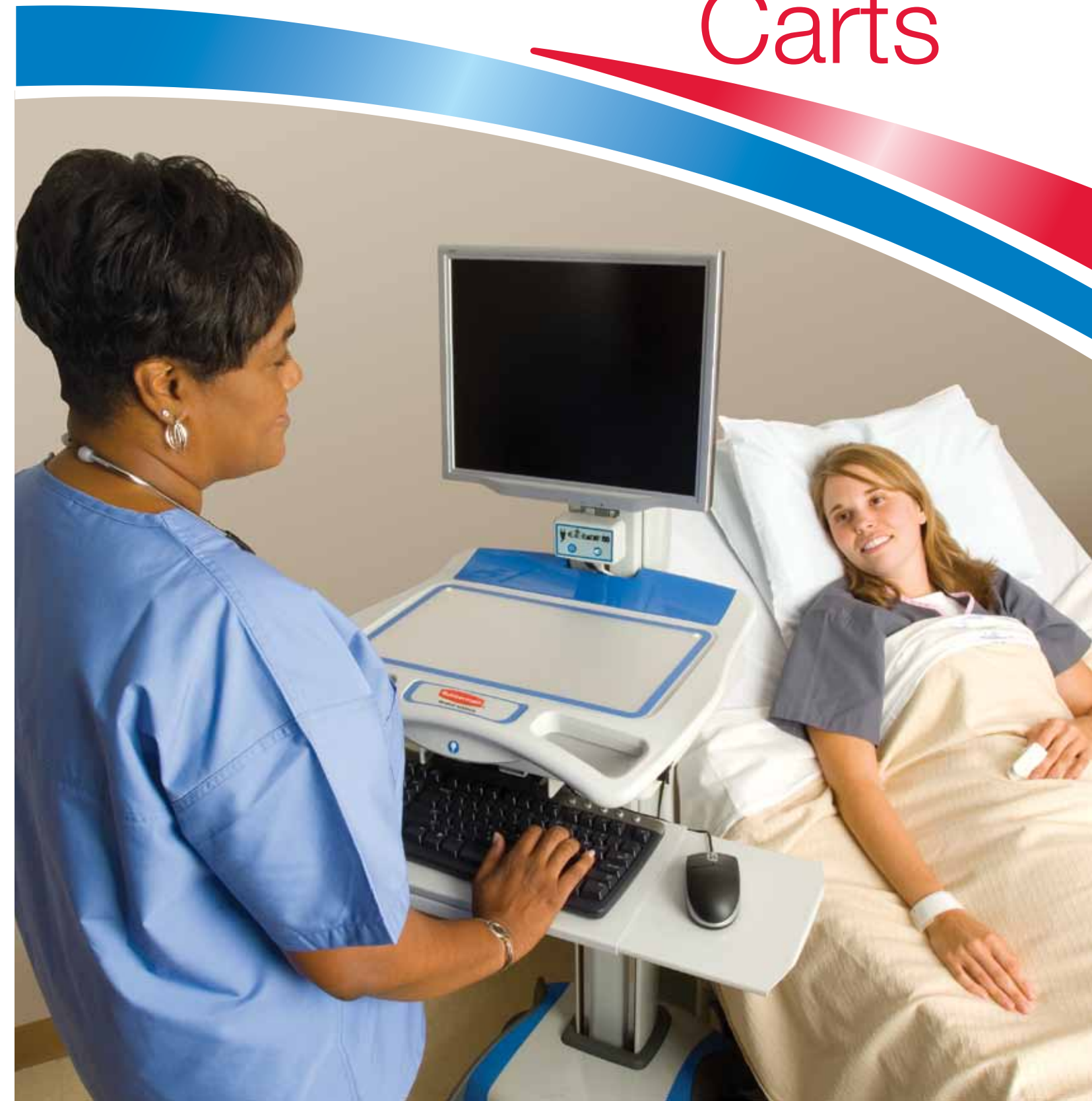



Medical Solutions

NEW! Lithium Batteries

- Longer Battery Life
- Shorter Charge Times
- Lighter Weight

Clinical Presence

Rubbermaid Carts are designed to promote a professional appearance specific to healthcare environments with internal cable management and concealed technology.

- Lightweight design, starting at 70 lbs
- Smooth, non-porous surfaces for easy cleaning and durability
- Open design prevents kicking the base
- Compact 17" x 17" base, contoured to 14.5" wide at center for easy maneuvering and storage
- Lithium Iron Phosphate (LiFePO4) Battery option reduces weight for increased maneuverability



Cart below shown with Medication Expansion Pack



Usability

Features designed for maximum maneuverability, smooth surfaces, and dedicated comfort grips improve clinician interaction with the product.

- Generous 21" w x 16" d work surface with raised edge to retain spills
- Document cover protects reference documents placed on the work surface
- Integrated keyboard light with auto shut-off feature for night time use
- Contoured handles on the front and back provide comfortable maneuverability from any position
- Accessories and storage options include but are not limited to: Wire Basket, Side Bins, IV Poles Vital Signs Mounts, Scanner/Printer Shelves, etc.

keyboard light



contoured handles



Workflow

Rubbermaid Medical carts are developed to improve workflow for nurses. These workflow features range from NEW lithium batteries to unique storage options.

- Medication Expansion Pack improves nurse workflow by providing organized, secure storage and capability to reconfigure drawers on the go
- New lithium battery options now available
 - Maximizes up-time
 - Minimizes charge time and overall weight
- 5 drawer sizes allow for 8 unique configurations to accommodate every storage need
- LCD Access Pack organizes essential supplies at the point of care
 - Medication Cup Dispenser
 - Dedicated tape storage
 - Alcohol wipe dispenser

Medication Expansion Pack



LCD Access Pack



Patient Interaction

Intuitive and simple adjustments over a wide range of motion allow clinicians to achieve eye contact with patients when using EMR software.

- Work surface with 16" of height range allow seated bedside interaction
- Premium monitor bracket allows monitor rotation for improved visibility while pushing cart and improved patient interaction at point of care.

keyboard tray



Patents Pending

Support Requirements

Easy to Select your Devices

Multiple power systems and battery options available; each with unique features tailored to fit a variety of workflow, technology and budgetary needs.


BATTERY OPTIONS	Rubbermaid Medical Solutions NON-POWERED	Rubbermaid Medical Solutions AC POWERED HIGH EFFICIENCY		Rubbermaid Medical Solutions DC POWERED HIGH EFFICIENCY
		LiFePO4	SLA	SLA
KEY BENEFITS	Lightest Weight Supports Self-Powered Devices	Longest Battery Life Span Shortest Charging Time Lightest Weight of Powered Carts Lowest Total Cost of Ownership	Low acquisition cost Upgradable to LiFePO4	Longest Run Time Upgradable to LiFePO4 Lighter Weight

Easy to Upgrade

Cart upgradeability extends the useful life of your carts.

POWER UPGRADE:

- ✓ Non-powered carts are pre-wired to accept power systems accommodating changes in technology and user preference
- ✓ SLA batteries are easily upgraded to Lithium batteries to minimize weight, support costs and maximize up-time

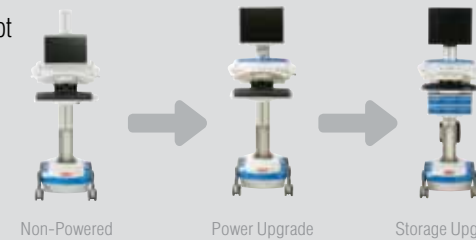


STORAGE UPGRADE:

- ✓ Pre-wired column allows for the addition of the Medication Expansion Pack supporting future bar-coding initiatives and changes in nursing workflow*

LIFE CYCLE MANAGEMENT:

- ✓ Maximize the return on your investment through the various phases of ownership
- ✓ Services include: pilot programs, customization, integration, battery management, preventative maintenance, and more

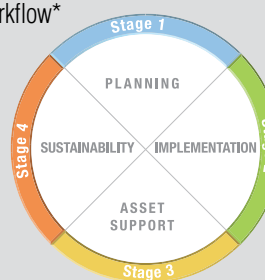

Easy to Maintain

POWER SYSTEM SOFTWARE

- ✓ Wirelessly monitors and logs system performance
- ✓ Email notification of low charge level prevents workflow interruption
- ✓ Proactive battery replacement notification

WARRANTY

- ✓ Next business day on-site service and parts as necessary
- ✓ Standard coverage:
 - 3 years mechanical
 - 2 years electrical
- ✓ Extended coverage available up to 5 years

Support Requirements