

M38 RX Mobile Medication Models



Self-Powered Device Medication Cart



Powered Medication Cart



Transfer and Exchange Carts

Mobile Medication Carts

NON-POWERED CARTS

- Configurations starting at 102 lbs
- Accommodates up to 17" notebook and tablet PC's
- External spiral cord recharges notebook
- Will accommodate LCD monitors if used in applications where cart is always plugged in

ALL POWERED CARTS

- Batteries: SLA AGM 35 or 55 AMP or Lithium Iron Phosphate (LiFePO4) 40 AMP*
- Certified to UL 60601-1 standard for patient point of contact
- Input 120V 60 Hz; Charger: 10 AMP SLA and 20 AMP LiFePO4
- Automatic transfer switch enables full use of PC without affecting battery charge time
- Power indicator at work-surface height with audible alarm
- Accepts LCD monitors up to 20"
- Auto shut down of PC prevents deep discharge
- Software tracks battery age and notifies when replacements are due
- Email setup sends cart specific critical alerts

AC POWERED CARTS

- Configurations starting at 124 lbs
- 3 AC power outlets
- External CPU configuration accommodates 2 ¼" to 4" wide CPU's

DC POWERED CARTS

- Configurations starting at 122 lbs
- High efficiency system provides more than 12 hours of run-time**
- 3 DC adapter power outlets
- Supplemental on-screen battery indicates time remaining

LITHIUM IRON PHOSPHATE (LiFePO4) BATTERY

- 40 AMP LiFePO4 provides 5 times longer life than 55 AMP SLA
- Charging Time: 2 Hours
- Weight: 14 lbs (65% lighter than SLA equivalent)

TRANSFER & EXCHANGE CARTS

- High capacity Transfer Cart (up to 144 drawers) for delivery from Pharmacy
- Unit based exchange cart with keyless lock system
- Removable Cassettes
- Flexible Drawer configurations
- Custom solutions available

*LiFePO4 for DC Carts available Q2 2011
 **Run times will vary depending on technology used and charging patterns

All M38 RX Models

GENERAL

- Automated Height Adjustment: 16" range of motion
- Contoured front and rear handle
- Primary Work Surface size is: 21" w x 16" d
- Primary Work Surface document protector: 16" w x 9.5" d
- Secondary slide out Work Surface size is: 8 ½" x 11"
- LED Keyboard light with auto shut-off
- Accommodates 8 3/4" d x 18" w keyboards
- Keyboard tray has left/right mouse pads, gel wrist rest, and easy access mouse holder
- Bin Accommodates barcode scanner heads up to 4x3"
- Casters: 4" precision bearing casters/2 locking
- T-track provides universal mounting point for accessories
- Smooth, cleanable surfaces (no ribbing)
- 2.5 ft hospital grade spiral cord extends to 8 ft
- Internal cable management
- Internal CPU storage: 15 ¼" w x 13 ¼" d x 3 ¼" h

MEDICATION ADMINISTRATION

- Electronic Keyless Entry
- PIN Code System: 0-9 digit, 128 codes per cart
- Network software for remote group PIN code update and audit
- Drawer Lock Time-out (seconds): 1-255 (default=120)
- Accommodates most bar code scanners

STORAGE

- 2 general purpose storage bins
- Storage for cables under work surface
- Wire basket (optional)
- Open side bin (2) (optional)
- Locking side bin (2) (Optional)
- Inner Drawer Dimensions (Maximum number per Cart)

- | | |
|---|---|
| ○ Small: (12) 4" w x 10 ¾" d x 2" h | ○ Medium: (8) 6" w x 10 ¾" d x 2" h |
| ○ Large: (4) 12" w x 10 ¾" d x 2" h | ○ Medium Double Deep: (4) 6" w x 10 ¾" d x 5" h |
| ○ Large Double Deep: (2) 12" w x 10 ¾" d x 5" h | ○ Envelope File Drawer: (1) fits 4-8 envelopes up to 11 ¾" x 9 ¾" |

SUITED FOR

- Initiatives
- Bedside documentation
- eMAR / Bar Code Medication Administration
- Units/Functions
- Med-Surg, Respiratory, ICU
- Centralized Medication Delivery (Cart-fill)
- Decentralized Medication Delivery (cabinet-based)
- Envelope or Unit Dose Delivery



SERVICE GUARANTEE

We are committed to providing best in class service to maximize your uptime. Our service guarantee includes:

- 2 hour response time
- Parts shipped next business day
- Help desk support
- Dedicated Technical Account Manager



Visit Rubbermaidmedical.com or call 888-859-8294 to see our complete line of medical products



Mobile Medication Carts





Medical Solutions

NEW! Lithium Batteries

- Longer Battery Life
- Shorter Charge Times
- Lighter Weight

Clinical Presence

Rubbermaid Carts are designed to promote a professional appearance specific to health care environments with internal cable management and concealed technology.


- Lightweight design, starting at 102 lbs
- Smooth surfaces for easy cleaning and durability
- Open design prevents kicking the base
- Compact 17" x 17" base, contoured to 14.5" wide at center for easy maneuvering and storage
- Lithium Iron Phosphate (LiFePO4) Battery option reduces weight for increased maneuverability.



slide out work surface



keyboard light

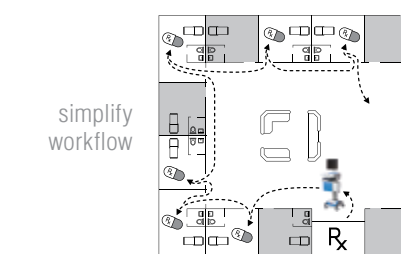


simplify workflow



LCD Access Pack (optional)



locking side bin (optional)

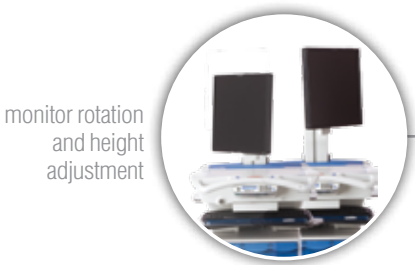

monitor rotation and height adjustment

Usability

Features designed for maximum maneuverability, smooth surfaces, and dedicated comfort grips improve clinician interaction with the product.

- Generous 21" w x 16" d work surface with raised edge to retain spills
- Secondary work surface slides out left or right
- Integrated keyboard and drawer light with auto shut-off feature for night time use
- Highly visible battery status indicator at work surface height
- Contoured handles on the front and back provide maneuverability from any position

Workflow

Nursing workflow is improved by providing organized, secure storage and the ability to reconfigure drawer on the go.

- 5 drawer sizes enable over 500 unique storage configurations to accommodate every storage need
- New lithium battery options now available
 - Maximizes up-time
 - Minimizes charge time and overall weight
- LCD Access Pack organizes essential supplies at the point of care
 - Medication Cup Dispenser
 - Dedication tape storage
 - Alcohol wipe dispenser
- Locking Side Bin provides additional locked storage that is removable for supply retrieval in the Med Room

Patient Interaction

Intuitive and simple adjustments over a wide range of motion allow clinicians to achieve eye contact with patients when using EMR software.

- Electronic lift with 16" of height range allow seated bedside interaction
- Premium monitor bracket allows monitor rotation and height adjustment for improved visibility while pushing cart and enhanced interaction at the point of care
- Improve nurse workflow by providing organized, secure storage and capability to reconfigure drawers on the go

Patents Pending

Support Requirements

Easy to Select your Devices

Multiple power systems available; each with unique features tailored to fit a variety of workflow, technology and budgetary needs.

BATTERY OPTIONS	NON-POWERED	AC POWERED HIGH EFFICIENCY		DC POWERED HIGH EFFICIENCY
		LiFePO4	SLA	SLA
KEY BENEFITS	Lightest Weight Supports Self-Powered Devices	Longest Battery Life Span Shortest Charging Time Lightest Weight of Powered Carts Lowest Total Cost of Ownership	Low acquisition cost Upgradable to LiFePO4	Longest Run Time Upgradable to LiFePO4 Lighter Weight

SLA batteries are easily upgraded to Lithium batteries to minimize weight, support costs and maximize up-time

Easy to Integrate technology

CABLE KEEPER:

- ✓ Dedicated cable storage for USB hubs, power bricks, and excess cables. All cables are concealed for a clean outward appearance.



CABLE CLIPS:

- ✓ Dedicated removable cable clips conceal excess cables



Easy to Maintain

- ✓ Optional PIN Code Management Software provides simultaneous PIN code updates to all carts from a single point



POWER SYSTEM SOFTWARE

- ✓ Wirelessly monitors and logs system performance
- ✓ Email notification of low charge level prevents workflow interruption
- ✓ Proactive battery replacement notification

WARRANTY

- ✓ Next business day on-site service and parts as necessary
- ✓ Standard coverage:
 - 3 years mechanical
 - 2 years electrical
- ✓ Extended coverage available up to 5 years

Support Requirements